



PEPLYN AIR Filters

- air / gas filters
- meltblown polypropylene

PEPLYN AIR filter cartridges have been specifically designed to guarantee removal of particulate from gas streams.

They can be used to protect sterilizing grade filters in pressurized systems or in exhaust gas vent applications.

PEPLYN AIR is particularly suitable for:

- Inlet gas in the fermentation industry as protection to sterilizing grade filters where polypropylene media is preferred
- As protection to sterilizing grade filters in exhaust gas systems
- Vent applications
- Systems where high particulate loading is expected

PEPLYN AIR has the ability to be steam sterilized and has a broad range of chemical compatibility


Features and Benefits

- Cost-effective prefiltration
- Steam sterilizable
- Absolute micron rating range from 1.0 - 25 micron
- Graded density for excellent particle retention
- High flow rates and long life
- No release of particles even during system pressure fluctuations




Note: PEPLYN is a registered trademark of Parker Hannifin Corporation.

Performance Characteristics

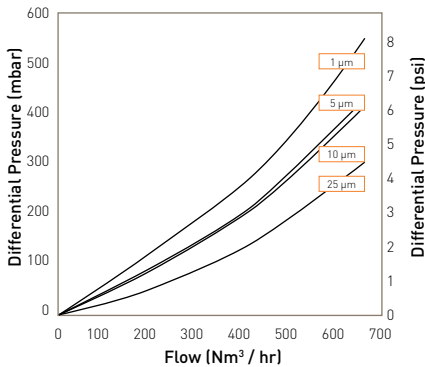

Flow rates for other sizes available upon request

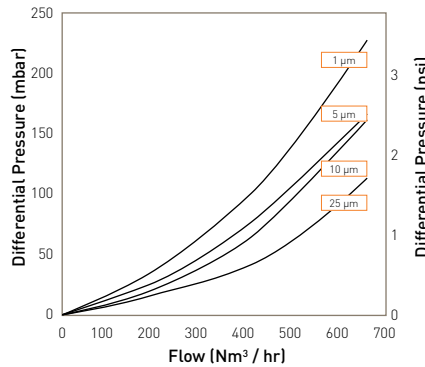

Flow rates for other sizes available upon request

**Cartridge flow rates @ 0 barg
10" Size (250 mm)**

**Cartridge flow rates @ 2 barg
10" Size (250 mm)**

Specifications

Materials of Construction

- Filtration Media: Meltblown Polypropylene
- Upstream Support: Polypropylene
- Downstream Support: Polypropylene
- Inner Support Core: 316L Stainless Steel
- Outer Protection Cage: Polypropylene
- End Caps: Polypropylene
- Standard o-rings/gaskets: Silicone

Biological Safety

Materials conform to the relevant requirements of 21CFR Part 177 and current USP Plastics Class VI - 121 °C and ISO10993 equivalents.

Recommended Operating Conditions

The maximum differential pressure in direction of flow (outside to in) is 3.5 barg (50.76 psig) at 20 °C (68 °F).

The maximum recommended continuous operating temperature is 50 °C (122 °F).

Effective Filtration Area (EFA)*

10" (250 mm) 0.49 m² (5.27 ft²)

*Varies with micron rating



Cleaning and Sterilization

PEPLYN AIR cartridges can be repeatedly in situ steam sterilized or autoclaved up to 142 °C (287.6 °F).

Determination of Micron Ratings

Particle removal efficiencies of PEPLYN AIR cartridges have been determined independently by challenging with a cut silica test dust, generated by BUS1701 dust injector used in conjunction with laser particle counters.

Micron Efficiency Ratings

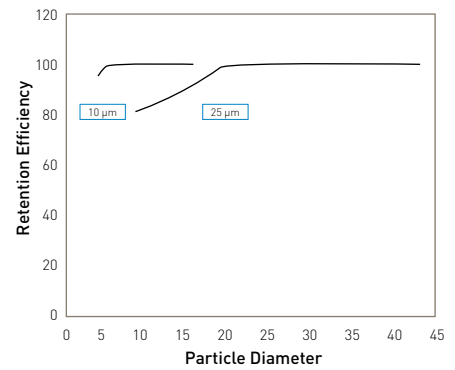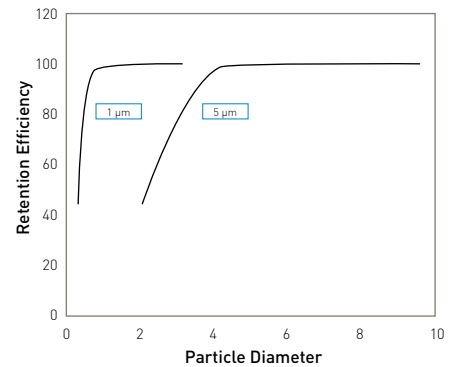

Ordering Information

ZCPH [] - [] - [] - []

Code Length (Nominal)	Code Micron	Code Endcap (10")	Code O-rings
B* 2.5" (65 mm)	1.0 1.0 µm	C BF / 226 Bayonet	E EPDM
A* 5" (125 mm)	005 5.0 µm	H UF Retrofit	P PTFE Encapsulated Silicone
K 5" (125 mm)	010 10.0 µm		S* Silicone
1 10" (250 mm)	025 25.0 µm	Code Endcap (Demi)	V Viton
2 20" (500 mm)		T TRUESEAL	
3 30" (750 mm)		Y Demi Stub	
		Z Demi A & B Std	

* Silicone o-ring supplied as standard without having to specify the 'S' code.

* Supplied in Packs of 3.

Europe: +44 (0)191 4105121 dhprocess@parker.com - North America: +1 608 824 0500 dhpsales.na@parker.com www.parker.com/dpharma

Parker dominick hunter has a continuous policy of product development and although the Company reserves the right to change specifications, it attempts to keep customers informed of any alterations. This publication is for general information only and customers are requested to contact our Process Filtration Sales Department for detailed information and advice on a products suitability for specific applications. All products are sold subject to the company's Standard conditions of sale.